


BUILT-IN PIEZO ULTRASONIC SCALER

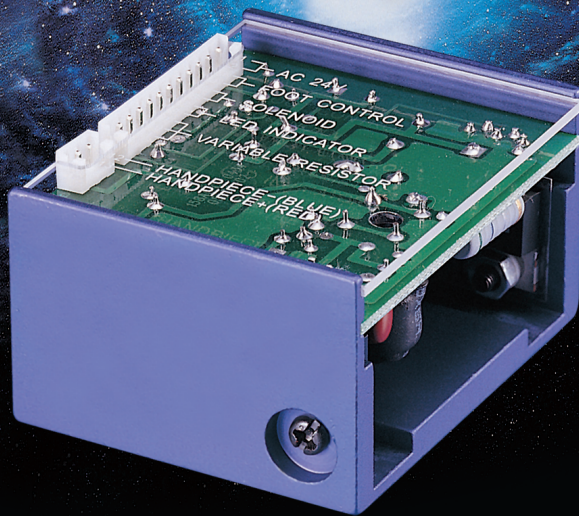

BUILT-IN PIEZO

- State-of-the-art piezoelectric technology is self-tuning and load-adjusting for smooth and effortless operation.
- Special TURBO function helps remove especially hard calculi by increasing power relative to the hardness.
- Streamlined compact size and slim lightweight handpiece.
- Innovative circuitry designed for greater stability during operation.
- Extremely gentle on the patient.
- Available selection of a wide variety of high quality tips to correspond with specific clinical applications.
- Unit is available for 100-240V 50/60Hz operation.
- Compatible with EMS type tips.

ISO 9001
ISO 13485



DENTAMERICA®

BUILT-IN

PIEZO ULTRASONIC SCALER

PIEZO SCALER TIPS



Reorder No. 88006
Perio



Reorder No. 88004
Universal B



Reorder No. 88005
Chisel



Reorder No. 88003
Universal A




Reorder No. 88002
Endo Chuck

SPECIFICATIONS :

- * INPUT: AC 24V 2A
- * FREQUENCY: 27,000 to 32,000 scaling strokes per second
- * DIMENSIONS: W65 x D50 x H33mm
- * WEIGHT: 330g
- * WEIGHT (HANDPIECE) : 40g


DENTAMERICA®

18688 E. San Jose Ave.
City of Industry, CA 91748 U.S.A.
Phone: (626) 912-1388
Fax: (626) 913-0510
info@dentamerica.com
www.dentamerica.com
www.dentamerica.asia

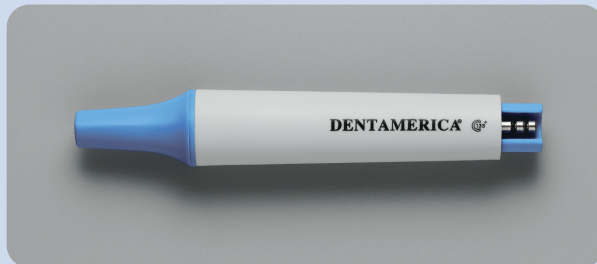

Reorder No. 88015
Handpiece

Autoclavable handpiece with piezo ceramic crystals and titanium tip socket structure.

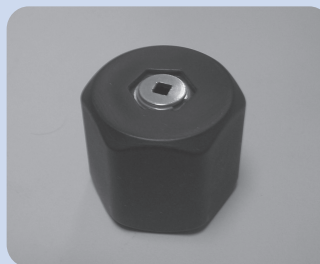

Reorder No. 88009
Torque-limiting
Tip Wrench

Special autoclavable wrench prevents damage to scaler tip when tightening.

CONTENTS

Control Box	x 1	Scaler Tip	x 3
Piezo Handpiece	x 1	Torque-limiting Tip Wrench	x 1
Handpiece Cable	x 1	Variable Resistor	x 1
11PIN Connecting Wire	x 1	Knob	x 1
2PIN Connecting Wire	x 1	Instructions	x 1

DISTRIBUTOR