

ENDOMAX[®]

CORDLESS ENDODONTIC HANDPIECE

ENDOMAX[®]



Order No.178



- Cordless convenience offers complete portability.
- Auto-stop and auto-reverse mechanisms for patient safety.
- Heavy-duty high torque motor.
- Reciprocating setting available without use of additional attachment.
- Ergonomic design provides comfort and easy operation.
- Features five torque settings.
- Automatic forward/reverse functions.
- Nine rotation speeds available for different applications.
- High capacity lithium-ion battery provides stable and prolonged usage.
- Offers quick and easy coupling.
- Low battery warning light.

ISO 13485




DENTAMERICA[®]

ENDOMAX[®]


CORDLESS ENDODONTIC HANDPIECE

Auto-Reverse Forward

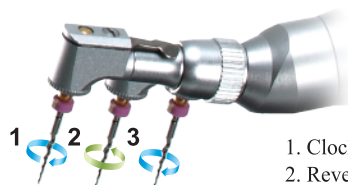


1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation continues forward.

Reciprocating Mode

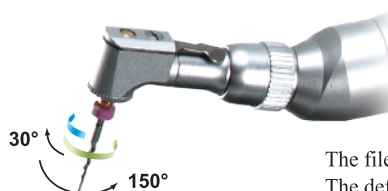

The file will reciprocate in a back and forth motion. The default rotation is 30° clockwise followed by 150° counterclockwise.

Auto-Reverse Stop

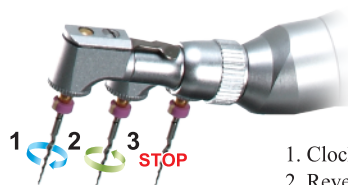

1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation stop.



Activation Button

Auto-Reverse Forward

Reciprocating Mode

Auto-Reverse Stop

Reverse Rotation Indicator

RPM Setting Button

Torque Setting Button

Low Battery Indicator

SPECIFICATIONS

INPUT	AC 100V~240V,50-60Hz
FUSE	1A 250V
BATTERY TYPE	3.6V 2040mAh
DIMENSIONS (HANDPIECE)	L192 xW26 xH29 mm
DIMENSIONS (CHARGING STATION)	L105xW85 xH59 mm
WEIGHT (HANDPIECE)	150g
WEIGHT (CHARGING STATION)	125g

DENTAMERICA[®]

18688 E. San Jose Ave.
City of Industry, CA 91748 U.S.A.
Phone: (626) 912-1388
Fax: (626) 913-0510
www.dentamerica.com
www.dentamerica.asia

DISTRIBUTOR

DENTAMERICA[®]